

GMFCS
Levels II–IV

kg
max. 65





41–98 cm

30–88 cm



I am your anterior walking aid, and will assist you in your very first attempts to walk.

I can hardly wait to get started with you! I support your balance and help train your torso muscles at the same time. You can even adjust my wheelbase to suit your needs: make it shorter and I become super manoeuvrable, so you can move me more easily. Choose a longer wheelbase and I give you more stability and support. And if you find it tricky to walk in a straight line, my direction locks are always ready to help.

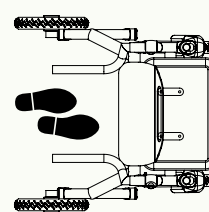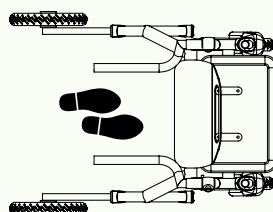



Handles

can be adapted in height and width.

Handbrake

with extra small brake levers with reduced reach.

Aluminium frame

highly compact; sizes 2+3 can be folded.

Regenerative brakes on the rear wheels

reduce the speed so that you can concentrate fully on walking.

Storage basket

provides space for toys or drinks.

Direction locks

support you in walking straight ahead.

PU tyres


are highly resistant.



Three dimensionally-adjustable forearm supports
stabilise your gait pattern.



Optional parking handbrake
provide you with more safety and control – also in combination with the forearm supports.

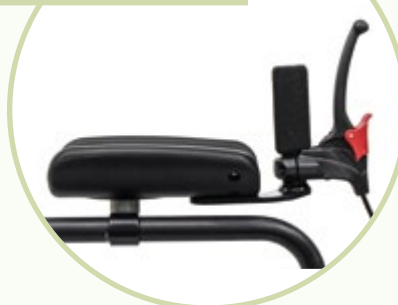

Frame colours
size 0+1



Frame colours
size 2



Frame colours
size 3





data	Size 0	Size 1	Size 2	Size 3
Handle height	30 - 47 cm	38 - 59 cm	53 - 75 cm	63 - 88 cm
Height of forearm supports	41 - 58 cm	49 - 70 cm	63 - 85 cm	73 - 98 cm
Handle spacing	18 - 32 cm	22 - 36 cm	30 - 44 cm	36 - 50 cm
Seat height	25 cm	30 cm	45 cm	
Overall width	44 cm	48 cm	62 cm	71 cm
Total length	36,5 cm (14,2 in)	42,5 cm	57,5 cm	61 cm
Overall length (telescopic wheelbase)	39 - 47 cm	45 - 55 cm	60 - 74 cm	65 - 79 cm
Max. load	20 kg	30 kg	50 kg	65 kg
Weight	4,2 kg	5 kg	7,4 kg	9,1 kg
HMV-No.	10.46.02.3049			



You can find more information on accessories in the price list.

