

Great sitting experience.

GMFCS
Level III-Vkg
max. 125

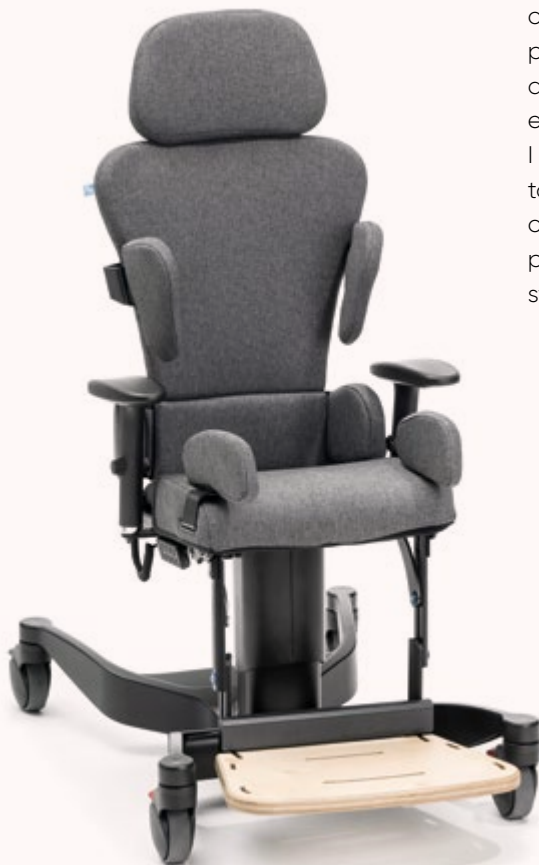
(-) 5° - 25°

(-) 20° - 20°

30 - 55 cm

40 - 55 cm

28 - 48 cm

**I am your grown up therapy chair.**

I'm your daily companion for schoolwork and comfort. I ensure you can sit with perfect posture no matter how tall you get. With my discreet design, I blend in seamlessly with everyday chairs, so you never feel out of place. I keep you positioned just right, allowing you to focus entirely on your tasks. My sturdy central column and extended wheelbase provide exceptional stability, making me strong enough to handle heavy use.

new



New tool-free adjustable version for headrest, chest pads and knee angle.



New size
We now offer **maxi.** in size 1b as well.



Front directional lock
are now included as standard and stabilise the straight-running performance.



Back unit

is stable and can withstand high loads.

Thigh guides

for optimum positioning.

Electric motor with hand control
for electrical adjustment of height and seat and backrest angle.

Castors

have a diameter of 100 mm and are equipped with two brakes at the rear.

Head rest
is upholstered and adjustable in height and angle.

Lateral pads
enable a stable and upright sitting position.

Arm rests
are width, height and depth adjustable and also can also be folded down.

Seat
is anatomically shaped for a comfortable sitting position.

Base frame
is compact and robustly manufactured.



High level of safety
thanks to a solid and stable base frame.



Hard-wearing cover
in modern colours with excellent quality.



The footrest
can be adjusted in height and knee angle and is permanently mounted to the seat unit.





new

data	Size 1	Size 1b	Size 2	Size 3
Seat depth	30 - 42 cm	38 - 48 cm		45 - 55 cm
Seat width	28 - 38 cm		36 - 48 cm	
Short back height	42 - 60 cm		48 - 60 cm	
Long back height	52 - 70 cm		58 - 70 cm	
Seat height	40 - 55 cm			
Lower leg length	27 - 45 cm		35 - 52 cm	
Total width	60,5 cm			
Length of base frame	84,5 cm			
Seat angle	(-)20° - 20°			
Back angle	(-)5° - 25°			
Max. load	125 kg			
Weight	45 kg			
HMV - No.	26.11.05.0051			

maxi. gives you greater independence with electric adjustment of the seat height, back angle and tilt-in-space functions.



You can find more information on fabrics
on pages 156 - 159.

You can find more information on accessories & harnesses
in the price List.

